



CAUTION

Before use, thoroughly read MAXAIR® CAPR® P/N 03521015 User Instruction (received with all CAPR Helmets and available at www.maxair-systems.com) for details.



REF 2099-05
(1 Set includes a right and a left Tab;
3sets/box)

Intended Use

The MAXAIR Liner Side Tabs are replacement side tabs for MAXAIR CAPR Helmet Liners.

Specifications

Material

White ABS

Recommended System Temperature Limits

Use/Handling: 0°C to 54°C at a maximum 80% Relative Humidity.

Charging: 0°C to 45°C at a maximum 80% Relative Humidity.

Storage: 0°C to 35°C at a maximum 80% Relative Humidity.

Regulatory

NIOSH

Symbol Definitions



Warning, Caution, or Note



Catalog Number



Part Number



WARNING

Failure to follow User Instructions P/N 03521015 and the instructions contained herein may be hazardous to the user's health.

Use only if package is received unopened and contents are undamaged. If damage is noted, contact the shipper for replacement or repair.

Prior to using any MAXAIR® System or component, be sure to be familiar with the system's NIOSH approved configuration.

DO NOT use if any component is damaged. If any components are damaged or contaminated and therefore unfit for safe and effective use, they should be replaced immediately.

Only trained and experienced personnel who have read and understand the User Instructions should use MAXAIR Products.

The institution using this product in any application is responsible for determining the appropriateness of this equipment relative to regulatory requirements. Bio-Medical Devices Intl, Inc. does not recommend the appropriate systems for a particular institution or facility.

Use only MAXAIR Systems/ NIOSH approved compatible components.

NOT for use in atmospheres immediately dangerous to life or health (IDLH), and atmospheres containing less than 19.5% oxygen, or more than 25% oxygen.


Flammability Level I: fabric may burn if exposed to open flame.

Follow current local regulations governing biohazard waste to safely dispose of single use MAXAIR Products.

If you need more information, contact your BMDI Sales Representative, or call BMDI customer service at 1-800-443-3842.

Liner Side Tabs

Front Button (Long Tab) Rear Button (Short Tab)

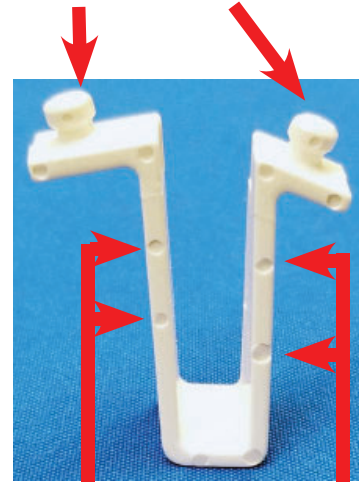

Circular hole indentations indicate this side to the outside of Liner.


Side Tab
(one on each side of Liner)

Two Liner slots on each side of Liner hold front and rear Buttons of Side Tabs



Liner Front and Rear Side Tab slots

Disassemble Liner Side Tabs

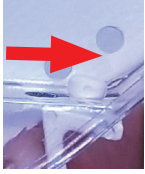

1. Push Rear Button backward and pull down to remove from rear slot.




2. Push Front Button inward and pull down to remove from front slot.



3. Side Tab removed from Liner. Repeat on other side.


Assemble Liner Side Tabs

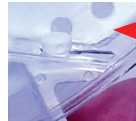

1. Insert Side Tab Rear Button into Rear Side Tab slot from underneath, then push all the way forward.



2. Insert Side Tab Front Button into Front Side Tab slot from inside, then push all the way forward (outward).

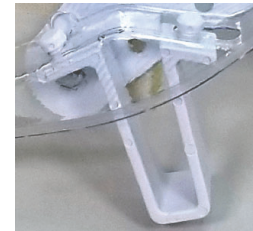

3. Properly installed Side Tab. Repeat for other side.